

Import sessions directly

aim


To import Viso sessions in The Observer XT. This makes use of the Noldus network communication protocol N-Linx.

prerequisites

- You have a license for The Observer XT 13 or a newer version.
- The computers with Viso and The Observer XT are in the same network and set up for communication with N-Linx. The procedure to set this up is present in the Viso Setup Manual.
- You recorded the sessions (Trainer), received session ownership or annotator rights or invitations to view the sessions or your user role is Administrator or Technician.

procedure


1. In The Observer XT, choose **File > Preferences**.
2. Open the tab **Viso settings**.

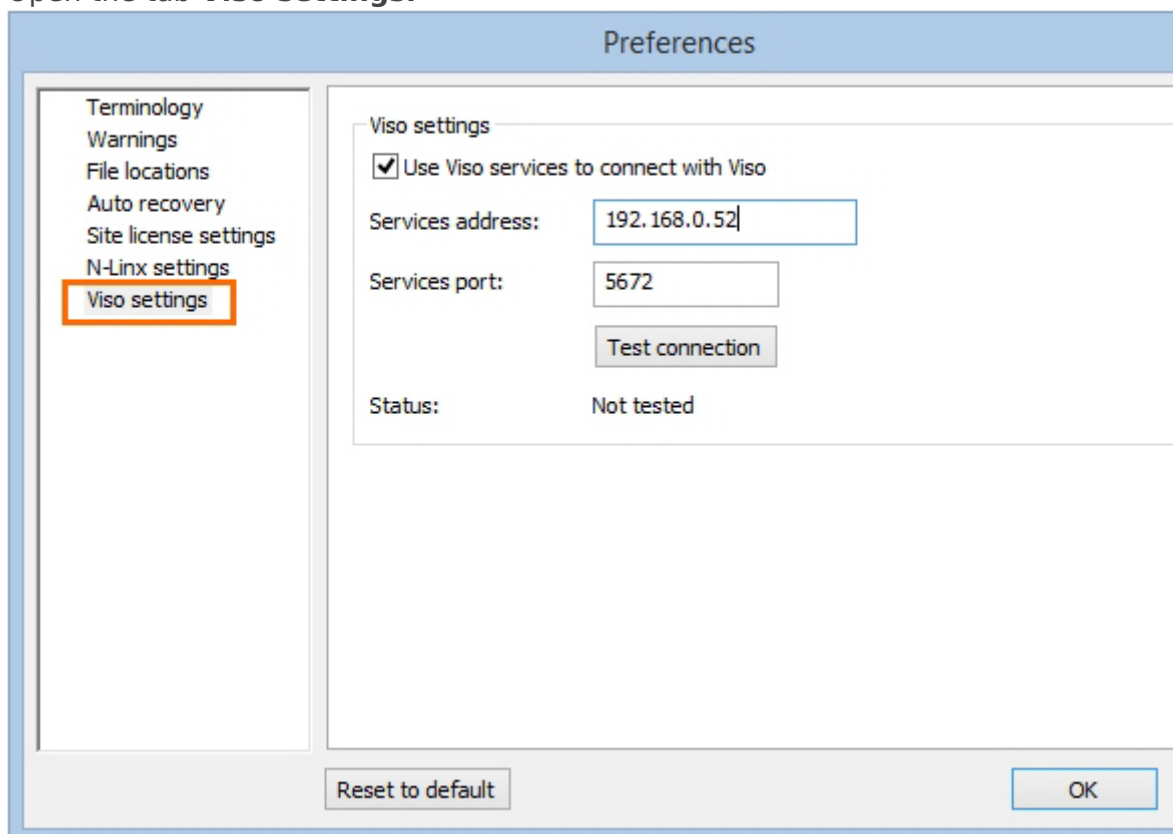

3. Select **Use Viso services to connect with Viso**.
4. In the **Services address** field, enter the computer name or IP address of the Viso Services computer. Ask your colleague who manages the Viso set-up for the details.
5. In the **Services port** field, leave the default port number 5672.

6. Click **Test connection**. If The Observer XT computer finds the Viso Services computer, the text **Connected** appears. Click **OK**.
7. Create a new project or open an existing one (**File> New/Open Project**).
8. Choose **File > Import > Viso Sessions**.
9. Enter your Viso login details.

NOTE If you are logged in in Viso, log out first. You can only import the sessions that you are allowed to view.

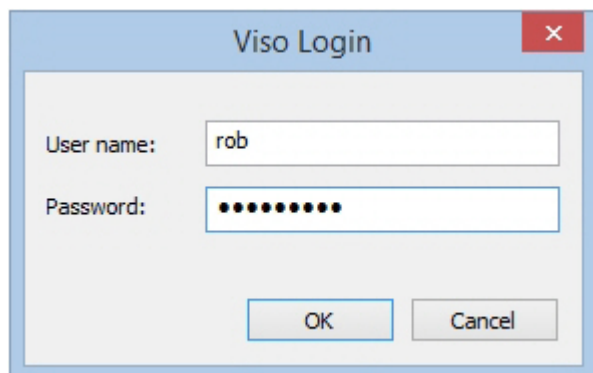

The image shows a 'Viso Login' dialog box with a title bar containing a close button (X). Inside the dialog, there are two input fields: 'User name:' with the text 'rob' and 'Password:' with a masked password represented by ten dots. At the bottom of the dialog are two buttons: 'OK' and 'Cancel'.

0. Select the sessions and click **Import**.

Import Viso Sessions					
	Name	Date	Duration	Location	Created by
<input checked="" type="checkbox"/>	Team skills training 1 debriefing	11/06/2017 12:00:00.29	30:00.18	Demonstration	Rose Bennett
<input type="checkbox"/>	Free play session	11/08/2017 11:19:00.08	39:25.31	Demonstration	Rose Bennett
<input type="checkbox"/>	Observation social interaction	11/08/2017 13:35:03.85	00:04.59	Demonstration	Rose Bennett
<input type="checkbox"/>	Mok counseling session 3	11/09/2017 09:30:00.09	29:59.99	Instruction roo	Rose Bennett
<input checked="" type="checkbox"/>	Team skills training 2 debriefing	11/08/2017 12:30:00.35	29:59.73	Instruction roo	Mr. R. Adams

NOTE Each session is imported as a separate observation. The videos are copied to the Video Files folder of The Observer XT. Markers are added to the coding scheme and remarks appear as comments in the Event log (see also [How Viso annotations appear in The Observer XT](#)). The observation is named after the imported session, with the location name and the date and time the session was created added to the filename.

NOTE If you have a license for The Observer XT 12.5 or earlier, you cannot import sessions directly into Viso. You must [Export sessions in Viso and import them in The Observer XT](#).

Revision #3

Created 15 January 2024 11:38:52

Updated 11 March 2024 14:26:38