

Record video


When you selected all cameras, you see a preview in the main MediaRecorder window.

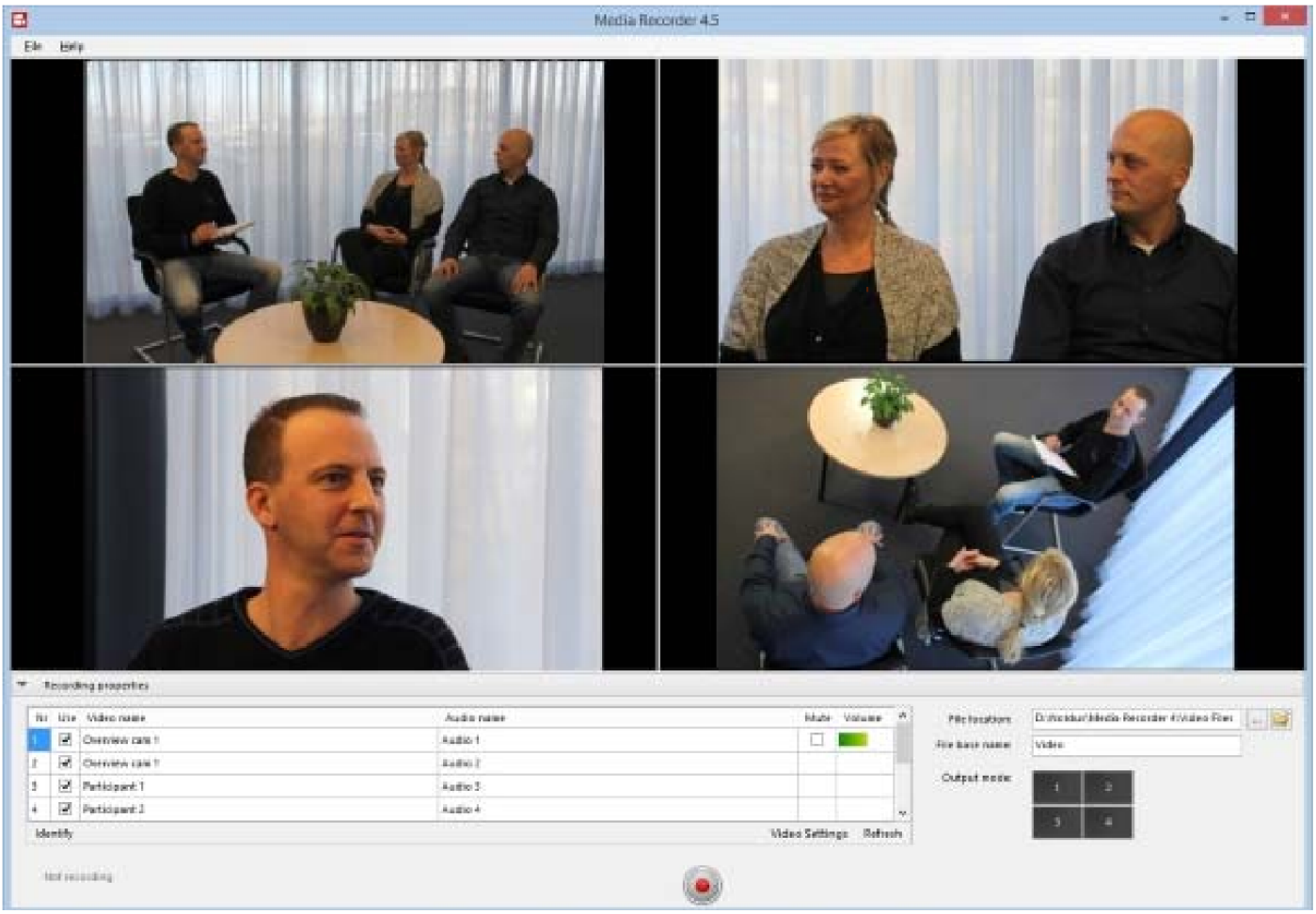


TIP If you selected multiple microphones, it can be useful to mute all but one audio channel in the preview. Select the audio source in the **Mute** column of the table under the preview. Muting audio in the preview does not influence the recording.

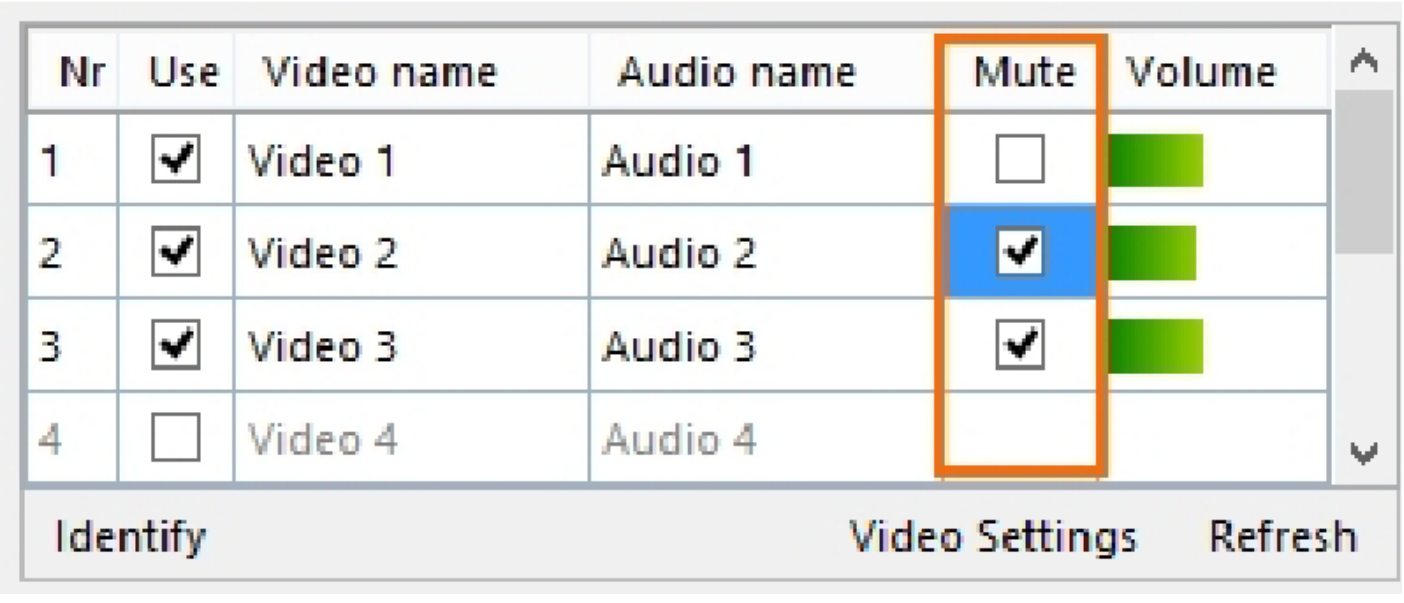

To start recording, press the **Start recording** button.

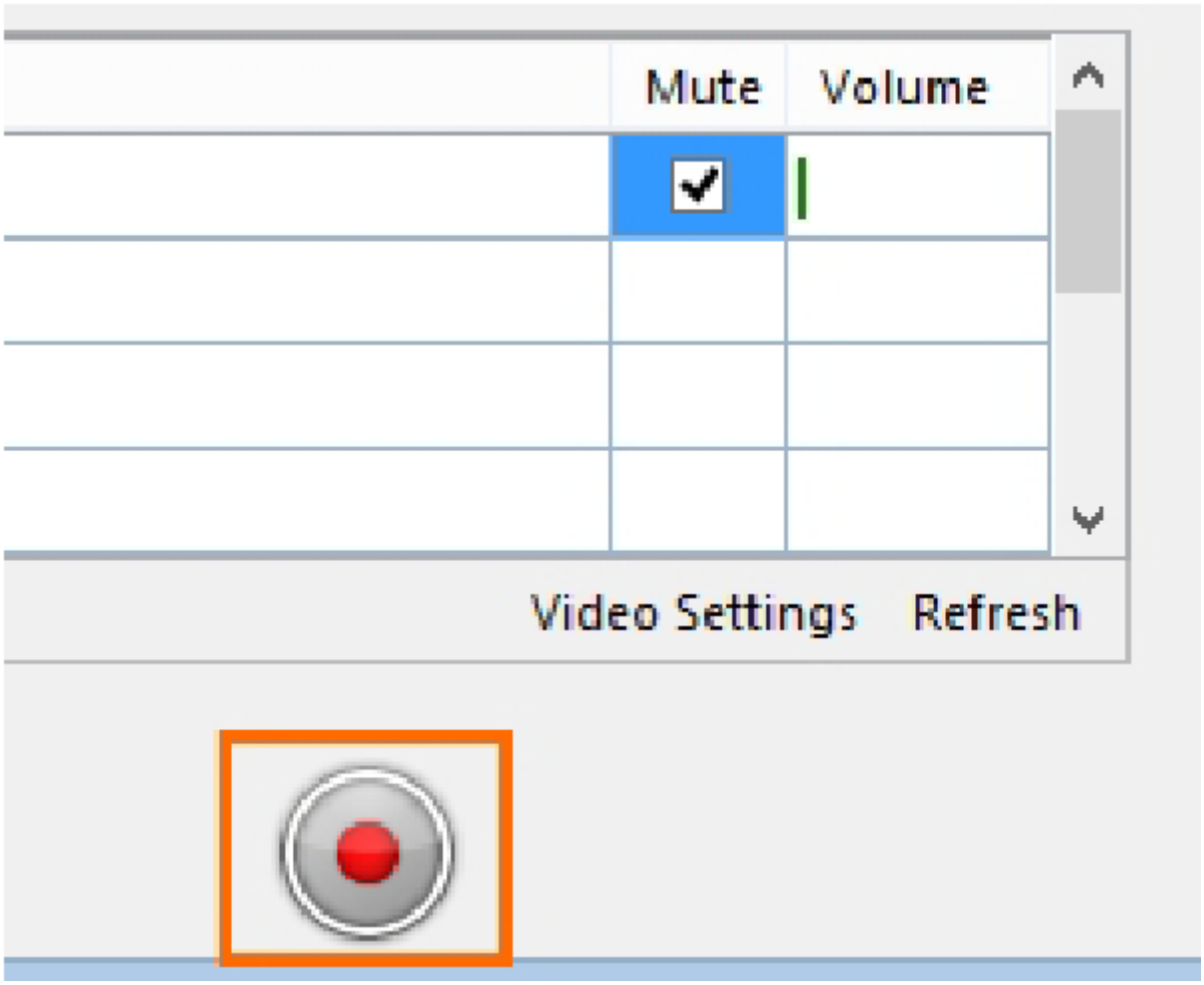

A pop-up window will appear, informing you about Noldus's Data Transparency Policy. Read it carefully before clicking the **I will inform the participants** button. You can also check the box **Don't show me again** to prevent further occurrence.

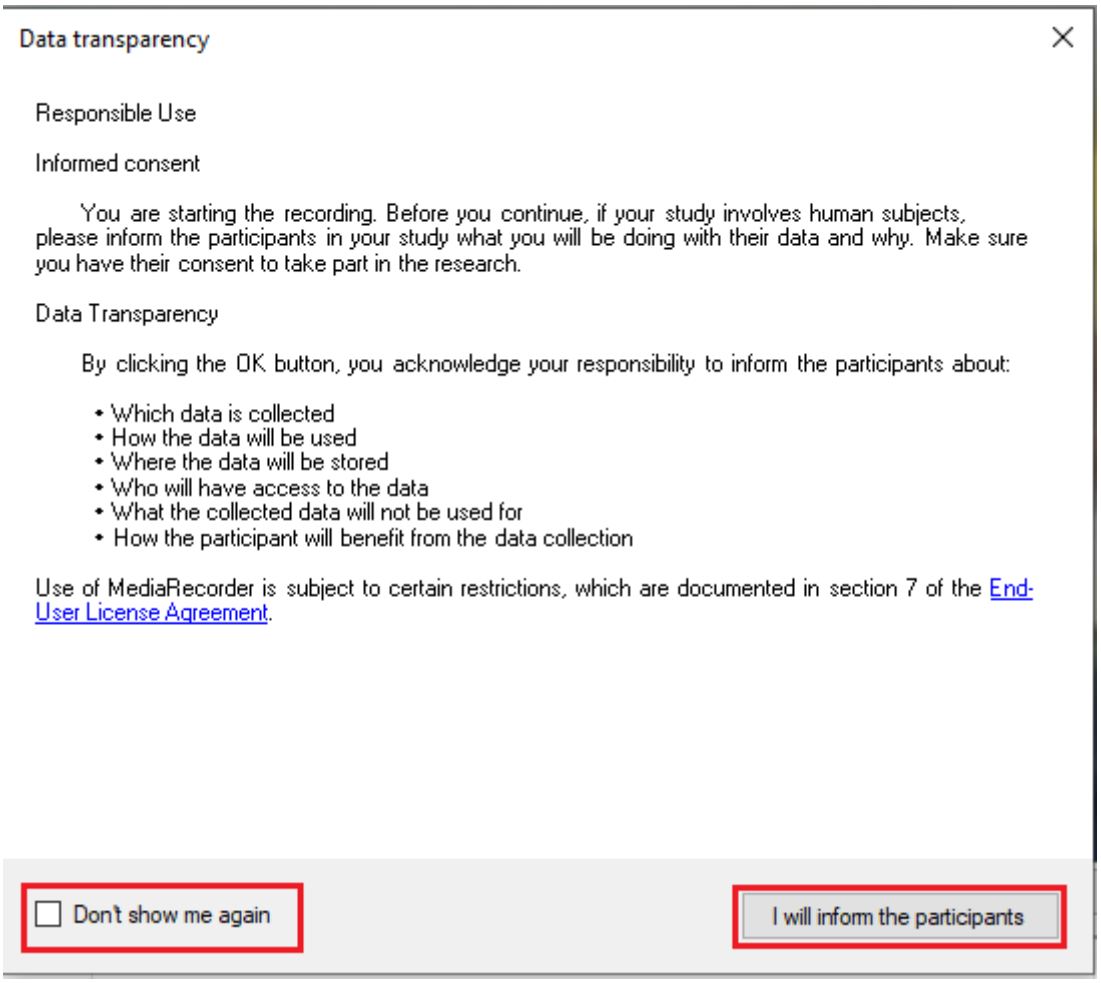


The button then changes into the **Stop Recording** button and the Recording timer starts running.



Press the button again to stop recording. Recording stops automatically when you specified a recording time .

how to access the recorded videos

To quickly access the recorded videos, press the **Open target folder** button in the main MediaRecorder window.

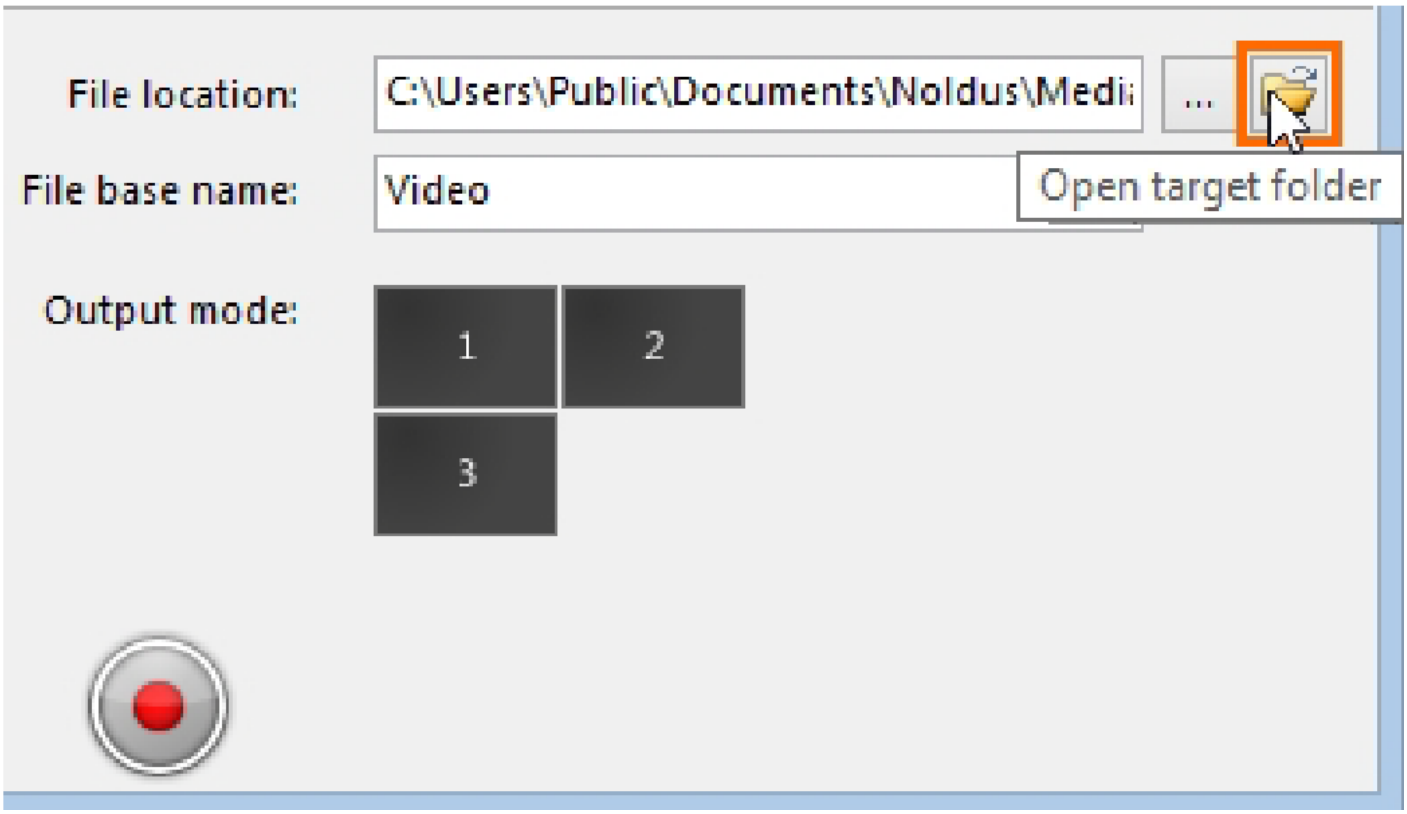

The video files are stored in the default video file location, or the location you specified. See To select the video file location.

Revision #8
Created 2023-11-29 14:49:07 UTC
Updated 2024-03-11 09:01:18 UTC