Create audio settings

Aim

To optimize the recording volume.

Prerequisites

- You have Axis IP cameras that are supported with MediaRecorder.

  See Supported IP cameras in Recording devices
- You followed the procedure in Ethernet card setup
- You connected all cameras.
- You gave the cameras an IP address.


  See Set camera IP address
- You gave the camera a password and selected the correct power line frequency.

  See Set camera password and select power line frequency
- You followed the procedure in Create camera settings

Procedure

1. Open the camera settings,

   See Access camera settings
2. Open the **Audio** tab.
3. Adjust the **Gain** slider and check in the **Input** section that the green field in the **Level** field is maximal with normal speech.

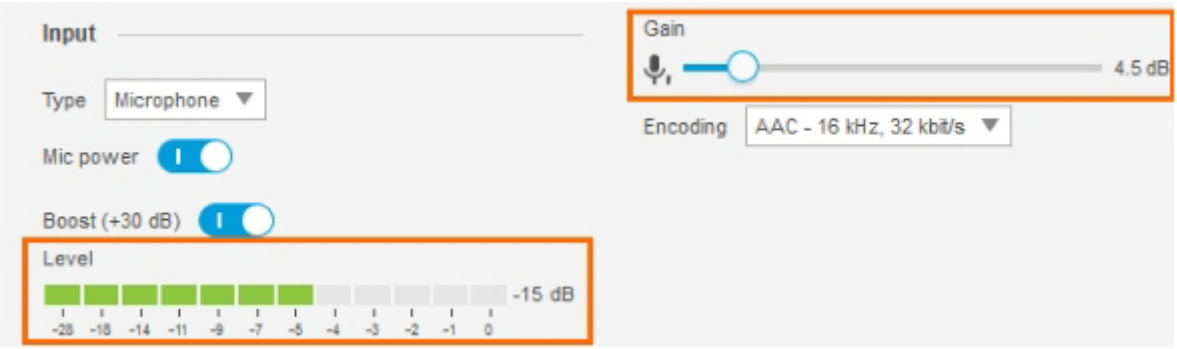

Note

You cannot record the audio of non-ONVIF devices in MediaRecorder. Optionally record audio from a separate microphone that is connected to the microphone input of the computer with MediaRecorder. The Axis cameras supported with MediaRecorder support ONVIF.

Revision #6
Created 2023-12-01 08:49:17 UTC
Updated 2023-12-20 12:47:08 UTC