# Set video file base name

## Aim

To select a word or phrase to start the video file name with.

The video file name consists of the **File base** name followed by:


- The video input name.
- The date and the time the videos were recorded.
- The number of the video camera.

Example – The **File base name** is Video, you changed the name of Video 2 to Webcam in the **Video Settings** window (**File** > **Video Settings**), the date is 11 May 2016, the time 10:54:43 and you create 2 videos simultaneously. The name of the second video will then be:

Video Webcam 5_11_2016 10_54_43 AM 2.mp4

## Procedure

Enter a name in the **File base name** field.

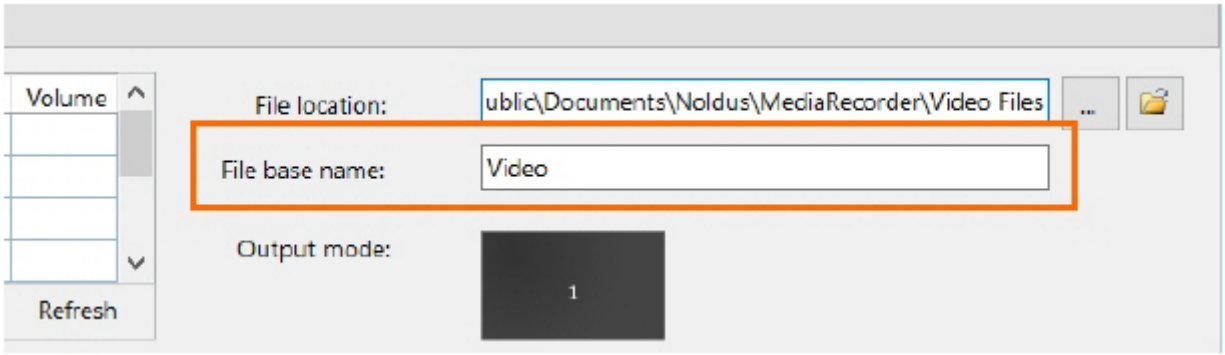


## Note

The date and time depend on your computer's regional settings. To change these, open the **Control panel** and then **Region and Language**. Select your preferences in one of the tabs.

---

Revision #3
Created 2023-11-29 21:37:59 UTC
Updated 2023-12-20 23:01:20 UTC