To get started


Main steps

To start working with MediaRecorder involves the following steps:

1. Connect all devices and, if applicable, switch them on.
2. If you have a hardware key - Insert it into a USB port of the computer.



3. Start MediaRecorder. The main window now opens. The first time it opens, the preview window is black. That is because you still have to select the cameras. Proceed with the next steps to do so.
4. Click Video Settings to select the cameras, microphones and, optionally, the output. See:
   Select a video device
   Select an audio device
   Video file options

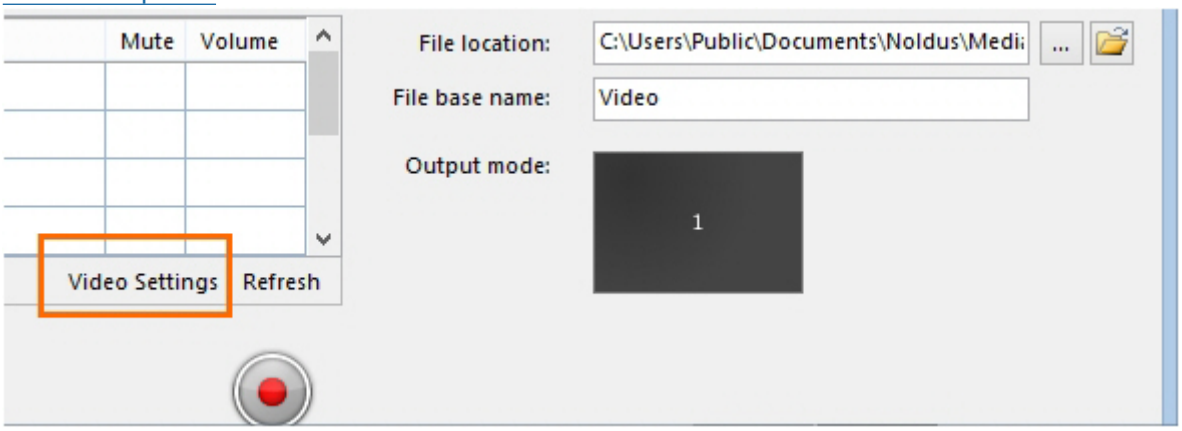

5. Record your videos

Revision #5
Created 2023-11-29 17:40:17 UTC
Updated 2023-12-20 12:47:08 UTC