Select multiple locations to record a composite session

aim


To select multiple locations for simultaneous recording. All recordings are combined in one session.

prerequisites

- The Technician account enabled multiple-location recording.
- Your user role allows session recording (Technician, Administrator, or Trainer).

procedure

1. In the **Locations** tab, select the locations.
2. Click the **Composite session** button on the toolbar.

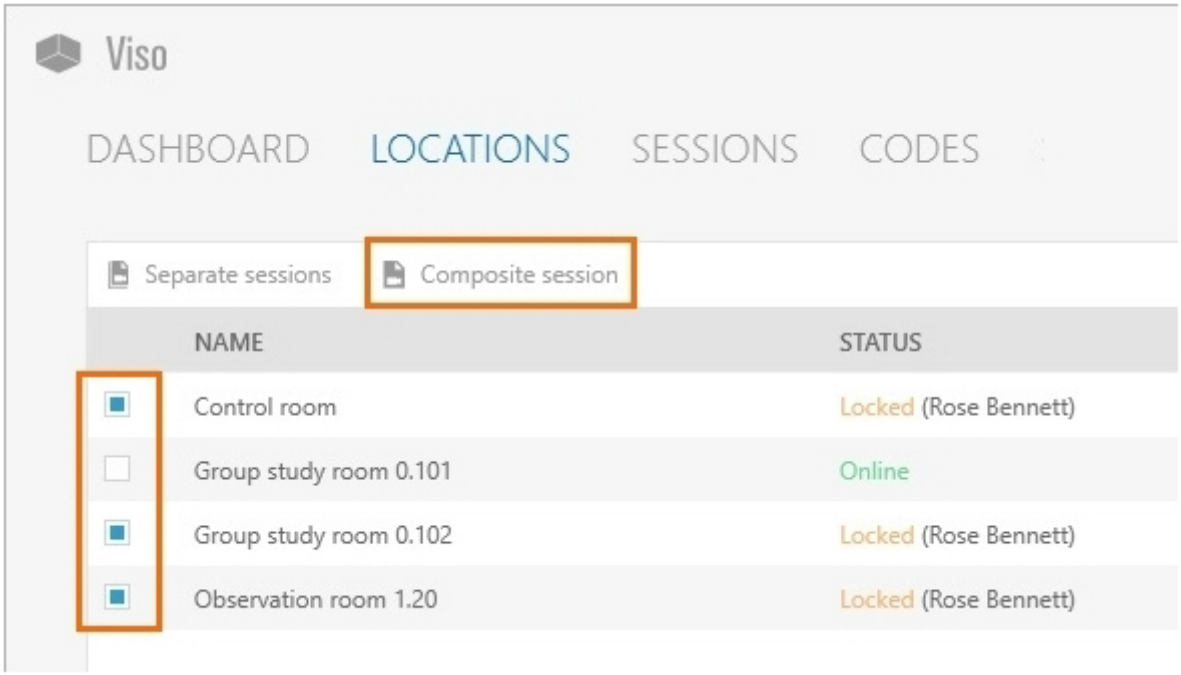

A preview of the cameras in the location that you selected first appears. The toolbar displays the carousel control.

See The carousel control and To open locations with the carousel control control to set up the locations

See the third chapter of this Reference (Prepare session recording) for the procedure. Follow the procedure for each location.

NOTE The Composite Session button is enabled once you selected multiple locations.

NOTE Composite session recording is supported with a maximum number of eight locations with two cameras each.

Revision #5
Created 2024-01-15 10:38:47 UTC
Updated 2024-05-31 11:56:53 UTC